

Sensitive adductor tendon

Example 1



Example 2



Example 3



With Your® Liner

Limb



3D model



Your® Liner



Indication


Patient complains about the compression applied by the liner on his medial adductor tendon



Typical issues

Off-the-shelf liner applies excessive pressure on adductor tendon

Off-the-shelf femoral liners are by essence cylindrical or conical. As off-the-shelf liners are designed to fit residual limbs of various lengths, they do not include a radius at the level of the adductor tendon, and thus apply extra pressure in this region, which can lead to discomfort, or skin irritations.

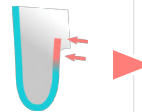

Key features



Anatomically tailored



Adapted compressions



Your® Liner includes radius at adductor tendon to relieve pressure

Your® Liner being designed for individual patients, a radius can be included in the 3D model to ensure the pressure applied on the adductor tendon is relieved, thus reducing the risk of the liner generating discomfort or skin irritations in this area.

How can Your® Liner help ?

